



Dedicated to Revolutionising the Power Chair Experience






QUANTUM®




W E L C O M E T O **QUANTUM**[®]

Quantum Rehab[®] was born out of the desire to delight customers with the most advanced, consumer-inspired complex rehab power wheelchairs and related technologies possible.

Clinically-focused and consumer-inspired, Quantum[®] products are engineered and manufactured to meet the widest range of complex rehab needs with functional, durable, stylised, high-performance designs to increase daily functionality and quality of life. At Quantum, consumer needs and wishes are the driving force. We're dedicated to not just meeting medical and clinical needs, but also quality-of-life needs. From the most advanced power seating for pressure management to USB ports, Bluetooth[®], and fender lights, no consumer need is overlooked.

Our Australian headquarters is based in Melbourne with additional warehouses in New South Wales, Queensland and regional Victoria. All our Quantum powerchairs are built in our Melbourne office, servicing all of Australia.





6	Edge[®]4
8	Outback[®] 4x4
10	4 Front[®] 2
12	Edge 3 Stretto[®]
14	Edge[®] 3
16	R-Trak[®]
18	Bases Continued
20	Bariatric Range
22	TRU-Balance[®] 4
24	iLevel[®]
26	TRU-Balance[®] 3
28	Seating & Positioning
30	Stealth Products[®]
34	Q-Logic 3
36	Accessories
38	Colour Guide
40	Specification Chart
43	PQ365 App
44	Service & Support
47	Stealth Resources

CONTENTS



TRU BALANCE[®] 4

POWER POSITIONING SYSTEMS



PROUDLY BUILT IN
AUSTRALIA



TRU-Balance[®] 3 & 4 Power Positioning Systems

Q-Logic 3 Advanced Drive Controls with Bluetooth[®]

OPTIONAL backup camera with 8.9 cm display can be mounted in multiple locations

OPTIONAL 12" power adjustable safe seat lift allows driving at speeds up to 7.2 km/h while elevated

4-pole, high torque motors

Stealth Products[®] cushion featuring exclusive Coolcore[®] technology

Front and rear full size LED lights for seeing and being seen (standard)

New battery tray to support 75 Ah battery

Comes standard with LED lighting package



"I LOVE HOW RESPONSIVE, AGILE, COMFORTABLE THE EDGE 4. IT'S AN AWESOME LOOKING MACHINE, I'M STOKED"

CAMERON IRWIN

EDGE® 4

MID-WHEEL DRIVE

The Edge® 4 is our most advanced model yet. With 180 mm (7") front and rear casters and powerful 4-pole motors, it delivers speeds up to 10 km/h and confidently tackles everyday obstacles. A new battery tray allows easier access and supports 75 Ah batteries for extended range. LED accent lights are integrated into the casters and fenders for greater visibility and style. Enhanced suspension and damping smooth out rough terrain, giving you comfort, control, and performance wherever you go.

SPECS

- **Maximum Speed:** 10 k/h (7.2 km/h when seat is elevated)
- **Weight Capacity:** 150 kg
- **Turning Radius:** 57.9 cm
- **Power Seat Elevation:** Optional 30.4 cm (12") iLevel

COLOURS

										
Pretty in Pink	Blue Jean	Evolution Orange	Lemon Crush	Bermuda Blue	Cosmic Purple	Forest Green	Ghost Gray	Ice Pink	Inferno Red	Royale Blue
										
Steel Blue	Night Train	Back in Black								



Automotive-grade SRS (Smooth Ride Suspension)





TRU-Balance[®] 3 Power Positioning Systems

Q-Logic 3
Advanced Drive Controls with Bluetooth[®]

OPTIONAL backup camera with 8.9 cm display can be mounted in multiple locations

OPTIONAL 10" power adjustable safe seat lift allows driving at speeds up to 5.6 km/h while elevated

Stealth Products[®] cushion featuring exclusive Coolcore[®] technology

4-pole, sealed, high torque motors


Front and rear full size LED lights for seeing and being seen (standard)

Off-road pneumatic 14" knobby tyres

Comes standard with LED lighting package

"IT'S GIVEN ME ACCESS TO THE BEACH AGAIN... I'M NOW ABLE TO GO OFF ROAD AND EXPLORE PLACES I HAVEN'T BEEN TO IN NEARLY 20 YEARS."

CAMERON IRWIN

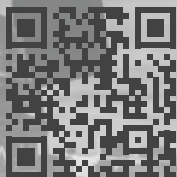

Scan to watch Cameron enjoy the Tasmanian coastline in his Outback 4x4 powerchair.

AVAILABLE WITH 10"
POWER SEAT ELEVATION

OUTBACK® 4X4

ALL-WHEEL DRIVE

The Quantum® Outback® 4x4 Powerchair is your ultimate adventure companion, ready to tackle almost any terrain with confidence. Designed for superior performance, this rugged powerchair features powerful 4-pole sealed motors and all-wheel drive for unbeatable traction and stability. Its robust construction ensures a smooth and comfortable ride on uneven surfaces, featuring advanced (SRS) suspension and off-road pneumatic 14" knobby tyres. Whether you're navigating gravel paths, bush tracks, or country properties, the Outback® 4x4 powerchair makes it a breeze.

SPECS

- Maximum Speed: 10 k/h (5.6 km/h when seat is elevated)
- Weight Capacity: 160 kg
- Turning Radius: 125 cm
- Power Seat Elevation: Optional 25.4 cm (10") elevation

COLOURS



Tan (Matte)

Sage (Matte)

Camo - Black

Camo - Desert

Camo - White



Automotive-grade SRS (Smooth Ride Suspension)





TRU-Balance[®] 4 Power Positioning Systems provide a choice of 10 or 30 degrees of power anterior tilt.

Q-Logic 3 Advanced Drive Controls with Bluetooth[®]

OPTIONAL backup camera with 8.9 cm display can be mounted in multiple locations

OPTIONAL 10" power adjustable safe seat lift allows driving at speeds up to 5.6 km/h while elevated

STC (Smart Traction Control)

Stealth Products[®] cushion featuring exclusive Coolcore[®] technology

Front and rear LED fender lights for seeing and being seen (standard)



Redesigned caster arms for improved durability

Automotive-grade SRS (Smooth Ride Suspension)

"RIDING IN THE 4FRONT 2 IS UNLIKE ANY CHAIR I HAVE EVER USED. THE 4FRONT 2 PUSHES THE LIMIT SET BY ITS PREDECESSOR BEYOND MY EXPECTATIONS."

CHRISTIAN BUDNEY



4FRONT[®] 2

FRONT-WHEEL DRIVE


FEATURING STC (SMART TRACTION CONTROL)

The 4Front[®] 2 is the next generation of the 4Front[®] Power Chair! New Smart Traction Control (STC), along with Smooth Ride Suspension (SRS), provide stability for smooth obstacle transitions and unmatched driving performance. Additional new features include 4-pole motors and redesigned caster arms for improved durability and aesthetics and reduced weight.

SPECS

- **Maximum Speed:** 10 k/h (5.6 km/h when seat is elevated)
- **Weight Capacity:** 136 kg
- **Turning Radius:** 65.8 cm
- **Power Seat Elevation:** Optional 25.4 cm (10") Safe Seat Elevation

COLOURS

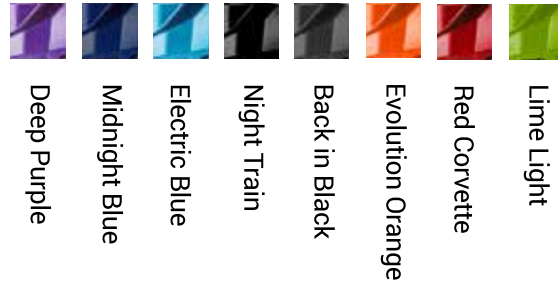

AVAILABLE WITH 10"
POWER SEAT ELEVATION



LED front and rear fender lights
(front: white, rear: red)



FUN FACT...

STRETTO IS THE ITALIAN WORD FOR NARROW



PROUDLY BUILT IN
AUSTRALIA



TRU BALANCE[®] 4

POWER POSITIONING SYSTEMS

TRU-Balance[®] 4 Power Positioning Systems provide a choice of 10 or 30 degrees of power anterior tilt.

OPTIONAL 12" iLevel[®] power seat elevation allows driving at speeds up to 5.6 km/h while elevated

Toolless removal of rear shroud with battery access door

Independent SRS (Smooth Ride Suspension)

OUR NARROWEST BASE IS ONLY
52 CM WIDE
WITH 12.5" DRIVE WHEELS

Q-Logic 3 Advanced Drive Controls with Bluetooth[®] and USB charger

Stealth Products[®] cushion featuring exclusive Coolcore[®] technology

Front and rear LED fender lights



Independent SRS (Smooth Ride Suspension)

Quantum's narrowest base



"TIGHT ROOMS, NARROW DOORWAYS, OUT-OF-REACH CABINETS AND VAN RAMPS...NOW I CAN GO WHERE I WANT TO GO, ALL BECAUSE OF THE STRETTO!"

BRYSON FOSTER

EDGE 3 STRETTO®

MID-WHEEL DRIVE

The Edge 3 Stretto® is the narrowest, most manoeuvrable complex rehab power base available in Australia! Aside from its narrow width, it is equipped with independent SRS (Smooth Ride Suspension), making it ideal for teenagers, and small adults. The Stretto is available with 31.8 cm (12.5") or 35.6 cm (14") drive wheels and optional 30.4 cm (12") of iLevel® power adjustable seat height.

SPECS

- **Maximum Speed:** 10 km/h (5.6 km/h when elevated with iLevel®)
- **Weight Capacity:** 136 kg (if increased weight capacity is required we recommend looking into the Edge4®)
- **Turning Radius:** 53.3 cm
- **Power Seat Elevation:** Optional 30.4 cm (12") iLevel

COLOURS



Deep Purple

Blue Jean

Back in Black

Green Machine

Red Corvette

Raspberry Beret

AVAILABLE WITH iLEVEL®
12"
POWER SEAT ELEVATION



Toolless access
to electronics





TRU-Balance[®] 4 Power Positioning Systems provide a choice of 10 or 30 degrees of power anterior tilt

OPTIONAL 12" iLevel[®] power seat elevation allows driving at speeds up to 7.2 km/h while elevated

Featuring Extra Stability Technology[®] for safer transfers, reach and other activities of daily living



Q-Logic 3 Advanced Drive Controls with Bluetooth[®]

Stealth Products[®] cushion featuring exclusive Coolcore[®] technology

Upgraded SRS (Smooth Ride Suspension)

USB charger (standard)

Front and rear LED fender lights for seeing and being seen (standard)



"THE EDGE 3 WITH ILEVEL HAS RESTORED MY INDEPENDENCE AT BOTH WORK AND HOME. ILEVEL HAS ALLOWED ME TO GET MORE INVOLVED AND START A DIALOGUE WITH PEOPLE, SOMETHING I WOULDN'T HAVE DONE BEFORE."

MERLISHA HENDERSON



Drive wheel color accents (standard)

EDGE[®] 3


MID-WHEEL DRIVE

The Edge[®] 3 offers upgraded SRS (Smooth Ride Suspension), which provides a smooth, comfortable ride and advanced stability. Optional iLevel[®] is available at 7.24 km/h with up to 30.4 cm (12") of power adjustable seat height.

SPECS

- **Maximum Speed:** 10 km/h (7.24 km/h when elevated with iLevel[®])
- **Weight Capacity:** 136 kg (if increased weight capacity is required we recommend looking into the Edge4[®])
- **Turning Radius:** 52 cm
- **Power Seat Elevation:** Optional 30.4 cm (12") iLevel

COLOURS

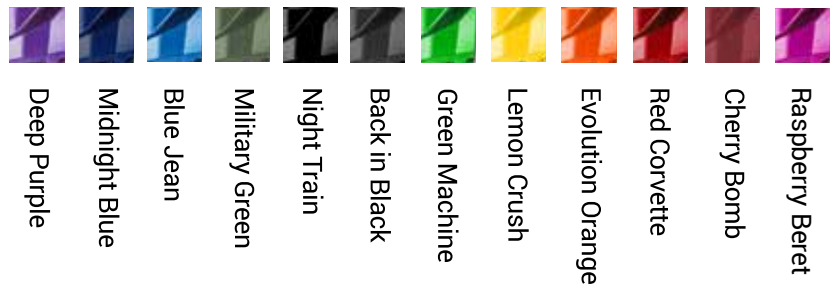

TRU-Balance[®] 3 Power
Positioning System

Q-Logic 3
Advanced Drive
Controls with
Bluetooth[®]

OPTIONAL 10" power
seat elevation allows for
driving at speeds up to
5.6 km/h while elevated

Featuring Extra Stability
Technology[®] for safer
transfers, reach and
other activities of daily
living

Stealth Products[®]
cushion featuring
exclusive Coolcore[®]
technology

Upgraded SRS
(Smooth Ride
Suspension)



Independent SRS (Smooth Ride
Suspension)

LED lighting in front caster
arms and rear fenders





R-TRAK

As a global innovator, we are creating a product that offers confidence and stability elevated at the speed of life. Our state-of-the-art facility has a robust team of experienced technicians who custom build each R-TRAK specifically to the individual consumer's needs.



Scan to learn about the design and development of the R-TRAK.

AVAILABLE WITH iLEVEL
 10"
 POWER SEAT ELEVATION

R-TRAK®

REAR-WHEEL DRIVE

The R-TRAK revolutionizes the rear-wheel drive power chair experience. Independent suspension on all four wheels combined with a compact base provide a smoother, more maneuverable driving experience. Standard features include fender lights, USB charger and drive wheel color accents. The R-TRAK is also available with our 25.4 cm (10") safe seat elevation, allowing the client to drive up to 5.6 km/h while elevated.

SPECS

- **Maximum Speed:** 10 km/h (5.6 km/h while elevated.)
- **Weight Capacity:** 150 kg
- **Turning Radius:** 65 cm
- **Power Seat Elevation:** Optional 25.4 cm (10") elevation

COLOURS



- Deep Purple
- Midnight Blue
- Electric Blue
- Night Train
- Back in Black
- Lemon Crush
- Red Corvette
- Lime Light




Optional Quantum backup camera can be mounted anywhere



Q6 EDGE® 2.0

MID - WHEEL DRIVE

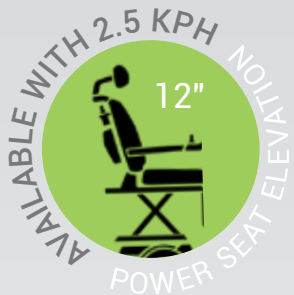


The Q6 Edge® 2.0 serves as Quantum's high-performance mid-wheel drive chair ideal as a dependable daily mobility companion. With standard high-speed (10 km/h) motors and ATX suspension, it delivers responsive control and stability—especially effective when custom spring packages are specified for lighter or paediatric clients. Optional iLevel® seat elevation (up to 30.5 cm / 12") at walking speed further enhances daily independence. The Edge 2.0 supports a wide range of seating and electronics options, offering exceptional value and versatility.


SPECS

- Maximum Speed: 10 km/h (2.5 km/h when elevated)
- Weight Capacity: 136 kg
- Turning Radius: 52 cm
- Power Seat Elevation: Optional 30.4 cm

COLOURS

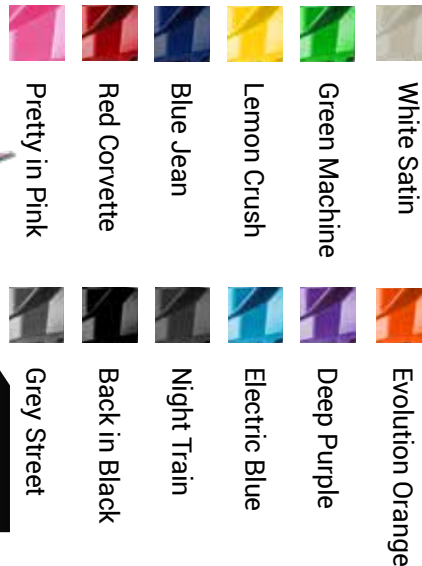


 PROUDLY BUILT IN AUSTRALIA

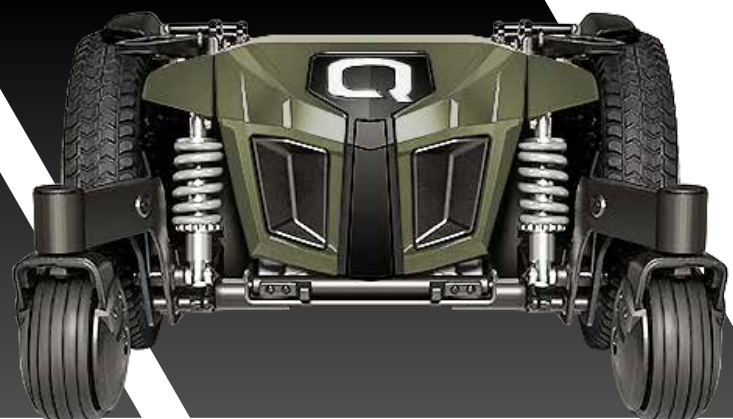

J623 2.0

MID - WHEEL DRIVE


The J623 2.0 is a high-performance mid-wheel drive chair engineered for active users seeking confident, responsive mobility. Powered by robust motors capable of speeds up to 10 km/h (6 mph), it offers smooth handling and strong outdoor performance. Standard TRU-Balance® 3 tilt (50°) provides comfort and positioning, while optional 25.4 cm (12") lift elevation supports greater access and independence. The VR2 controller enables two powered seating functions, offering simple and reliable controls for everyday driving and power positioning.

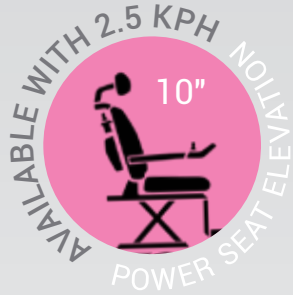

SPECS

- Maximum Speed: 10 km/h (2.5 km/h when elevated)
- Weight Capacity: 136 kg
- Turning Radius: 52 cm
- Power Seat Elevation: Optional 25.4 cm

COLOURS

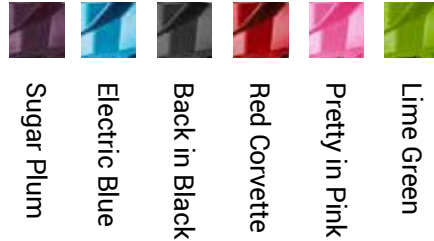

* ADDITIONAL J623 2.0 COLOURS INCLUDE: CARBON FIBRE, EMERALD GREEN, ORANGE CRUSH, SILVER BULLET, & YELLOW BLAZE

 PROUDLY BUILT IN
AUSTRALIA



Q6 EDGE® HD
204 kg Weight Capacity

QUANTUM® 1450
272 kg Weight Capacity



Q6 EDGE® HD

MID-WHEEL DRIVE

SPECS

- Maximum Speed: 7.6 km/h (5.6 km/h when elevated)
- Weight Capacity: 204 kg
- Turning Radius: 50 cm
- Power Seat Elevation: Optional 25.4 cm iLevel®



4 FRONT 2 HD

FRONT-WHEEL DRIVE

SPECS

- Maximum Speed: 8 km/h (4.3km/h when elevated)
- Weight Capacity: 204 kg
- Turning Radius: 65.7 cm
- Power Seat Elevation: Optional 25.4 cm

BARIATRIC RANGE

BUILT TOUGH

Quantum® bariatric power chairs all offer an innovative design with superior outdoor performance and superb indoor maneuverability. Our bariatric line of power chairs offers speeds as high as 8km/h and weight capacities up to 272kg.

All of our bariatric power chairs can be equipped with power positioning for greater comfort, such as TRU-Balance® HD Power Tilt. This seating system can assist with pressure relief. To provide additional support, Stealth® cushions and backs are available to provide lumbar support or coccyx and sacral relief. All Quantum bariatric power chairs are available with Q-Logic environmental controls featuring Bluetooth® technology.



QUANTUM® 1450

FRONT-WHEEL DRIVE

SPECS

- **Maximum Speed:** 7.2 km/h
- **Weight Capacity:** 272 kg
- **Turning Radius:** 67.3 cm
- **Power Seat Options:** Bariatric tilt

TRU-Balance® 4 Power Positioning Systems

WITH ANTERIOR TILT AND MEMORY SEATING

Maximum weight capacity of up to 150 kg with Edge® 4

Up to 8 customizable memory seating positions recall the saved positions of all actuators, offering the choice between sequential or synced movement for each seat function.

Seat elevation combined with anterior tilt provides optimal access.

Consumers may find it easier to perform ADLs such as: accessing sinks and cooktop surfaces, meal preparation and reaching items in cabinets.

TRU-Balance® 4 with Dual Actuator AFP shown on the Edge 3 Stretto®



Chest Bar




TRU-Balance® 4 Power Recline



TRU-Balance® 4 Power Posterior Tilt





DUAL ACTUATOR POWER ARTICULATING FOOT PLATFORM

Lowers to the ground to assist
with transfers

Articulation range of 0-7.25"

Elevation range of 85-165°

Only available on TRU-Balance® 4

TRU-BALANCE® 4

POWER POSITIONING SYSTEMS

FEATURES

- Anterior tilt options: 10 and 30 degrees
- Up to eight customizable memory seating positions
- Chest bar and knee blocks available (knee blocks and chest support required with 20 or 30 degrees of anterior tilt)
- Available on the Edge 3 Stretto®, Edge®4, Front® 2, and Edge® 3
- Optional Accu-Trac Advanced Tracking Technology
- WC19 crash-test approved occupied transit kit option
- Available seat sizes: 16-22" width and 16-22" depth

POWER RECLINE

85° - 165° of recline

POWER POSTERIOR TILT

0° - 50° of tilt

Maintains center of gravity shift

DRIVE SPEEDS IN ANTERIOR TILT

Up to 10° of anterior tilt, up to full lift: 3.2 km/h


10-20° of anterior tilt, up to 15.3 cm (6") of lift: 0.8 km/h

10-20° of anterior tilt, over 15.3 cm (6") of lift: drive lockout

Over 20° of anterior tilt: drive lockout



Knee Blocks



TRU-Balance® 4
with 12" iLevel®
Power Elevation



TRU
BALANCE® 4
POWER POSITIONING SYSTEMS



Safer Transfers



Improved Reach



Q
QUANTUM®

"COOKING IS SO MUCH EASIER NOW THAT I CAN REACH THINGS. I CAN'T IMAGINE MY LIFE WITHOUT ILEVEL."

STEPHANIE WOODWARD

iLEVEL®

POWER SEAT ELEVATION

Standing Height + Speed + Stability = Authentic iLevel (like no other)

iLevel® power seat elevation technology features a revolutionary new power chair design for a better experience. iLevel improves daily functionality and quality of life by allowing operation of the power chair with the seat elevated while at walking speed.

FEATURES

- 30.4 cm (12") of power adjustable seat height while at walking speed; maximum of 25.4 cm (10") height on the Q6 Edge® HD
- iLevel's patented Extra Stability Technology® enhances safety for transfers, reach, and numerous other activities of daily living
- Available on the Edge 3 Stretto®, Edge® 3, Edge® 4, Q6 Edge® 2.0, and Q6 Edge® HD
- Front and rear LED fender lights come standard with iLevel®

• • QUANTUM®
iLevel®
power chairs



Social interaction
at eye level

TRU-Balance® 3 Power Adjustable Seat Lift with iLevel® technology

Maximum weight capacity:
12" up to 150kg.
HD iLevel up to 181kg (Edge HD)

Provides extra stability to
ensure maximum safety
for transfers, reach, and
numerous other activities
of daily living

Maximum speed of 7.24
km/h while elevated at
30.4 cm of lift on the
Edge®3, Edge® 4

Seat-to-floor-height as low
as 44.45 cm on the Edge 3


LED fender lights
included



TRU-Balance® 3 Power Recline



TRU-Balance® 3 Power Tilt

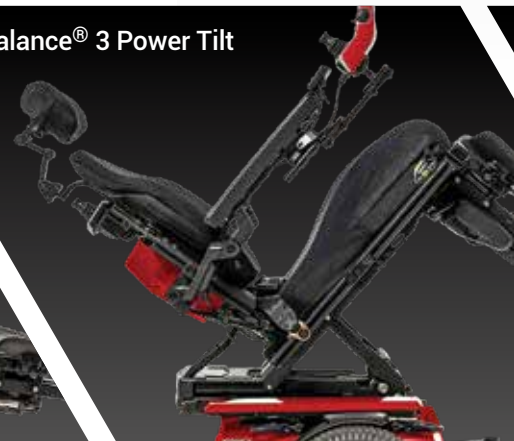

POWER ARTICULATING FOOT PLATFORM

Broad articulation range for enhanced adjustment

204 kg weight capacity

Unmatched adjustment, durability, and strength

Multiple calf pads, calf panel, and footplate options

TRU-BALANCE® 3

POWER POSITIONING SYSTEMS

FEATURES

- Fine adjustments can be made even while client is in the chair
- Seat width can be ordered independently of the back width at no additional charge
- Power tilt, power recline, and iLevel® are retrofittable

POWER RECLINE

85° - 175° of recline Available in matching back shroud colours

160 kg weight capacity

Available with power tilt and power adjustable seat lift

POWER TILT

50° of tilt Maintains center of gravity shift

Lower seat-to-floor height

Unmatched adjustment and durability

HEAVY DUTY POWER TILT & RECLINE

204 kg weight capacity Available in matching back shroud colours

Up to 50° of tilt 85° - 175° of recline

TRU-Balance® 3 Heavy Duty Power Tilt and Recline



TRUBALANCE® 3
POWER POSITIONING SYSTEMS



Synergy® Seating



Contoured Seating



MORGAN STEWARD ENJOYS THE SPEED, ELEVATION, AND RECLINE OF HIS ILEVEL® POWER CHAIR.

MORGAN STEWARD

SEATING

AND POSITIONING

SYNERGY® SEATING

Easily adjustable, this versatile design can grow with the user

Several width/depth adjustment packages ranging from 10" to 32"

Accepts a wide variety of after-market seats and backs

Heavy duty option provides a weight capacity of 272 kg and seating up to 32" x 28" (model specific)



Q6 Edge® 2.0 with Synergy® seating

BARIATRIC TILT

TRU-Balance HD Power Positioning Systems offers increased weight capacities and expanded seating sizes

Configurations with a 204kg weight capacity include HD power tilt and HD power recline

10" power adjustable seat lift is available with a 181kg weight capacity on the Q6 Edge® HD the 4Front® 2 HD.

Bariatric Tilt for configurations up to 272kgs is also available on the Q1450.

Quantum® paediatric seating with iLevel® 30.4 cm power seat elevation

Stealth Products® CUSHIONS


FEATURING EXCLUSIVE
coolcore®
CHEMICAL-FREE COOLING TECHNOLOGY



Convenient built-in
side/front pouches



Constructed with FRG
Sunmate flame-retardant
material

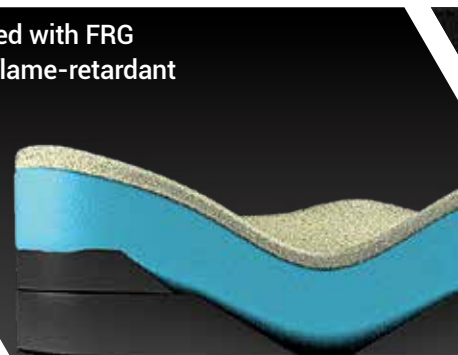

Q
QUANTUM®



30+ UPF



Dry



Chemical Free



Breathable Comfort

"While the other makers of cooling fabrics use a number of chemicals to enhance wicking, transportation and evaporation, Coolcore fabrics are chemical free. To make that possible, we changed the essence of cooling fabric technology right down to its very core."

COOLCORE®



STEALTH PRODUCTS®

CUSHIONS

Stealth Products® Cushions feature patented Coolcore® technology, a chemical free, four-way stretch material that offers moisture wicking, moisture transportation, and moisture evaporation properties. Our cushions also use a colour guide that identifies the suggested function family by a unique patch.

FEATURES





- Chemical free
- Regulates temperature
- Manages and distributes your body heat



Easy-pull rubber loop zipper handles





	TRU-Comfort 2 SPP	Essence SPP	Solution® SPP	Spectrum® Foam/Gel SPP
DESCRIPTION¹	 <p>Multi-layered foam Deep contour, CoolCore® cushion cover has “pockets” on sides</p>	 <p>Low profile Tru-Comfort 2, lateral and pre-ischial shelf Contoured, HD molded foam base, coccyx and sacral relief</p>	 <p>Standard viscoelastic foam top with twin cell aqueous viscous gel (heavier gel) insert in ischial well, mid-level contour</p>	 <p>SunMate FRG foam top (SunMate FRG XXSoft (blue); Medsoft (green)) with multi-layered foam insert or twin cell ultra-lightweight viscous gel insert in ischial well, sacral relief</p>
WEIGHT CAPACITY	14-15" W: 79.4 kg 16-21" W: 158.8 kg 22+" W: 204.1 kg	10-15" W: 79.4 kg 16-21" W: 158.8 kg 22" W: 204.1 kg	10-15" W: 79.4 kg 16-20" W: 158.8 kg 22+" W: 226.8 kg	14" W: 79.4 kg Foam & Gel 16-20" W: 158.8 kg Foam & Gel 22+" W: 294.8 kg Foam; Gel
OVERALL WEIGHT	1.8 kg	1.6 - 2.3 kg	2.4 kg	Foam: 1.4 kg Gel: 1.6 kg
HEIGHT	11.4 cm / 4.5"	8.9 cm / 3.5"	7.6 cm / 3"	7.6 cm / 3"
DIMENSIONS²	14x14" up to 24x24"	10x10" up to 22x24"	10x10" up to 24x24"	Foam: 14x18" up to 32x28" Gel: 14x18" up to 24x26"

ADI Q-Back 2.0



The ADI Q-Back 2.0 from Stealth Products® is available on TRU-Balance® 3 and TRU-Balance® 4 Seating Systems and is offered in three configurations:

- Active Contour ADI Q-Back 2.0
- Deep Contour ADI Q-Back 2.0
- Extra Deep Contour ADI Q-Back 2.0
- Coolcore® cover comes standard

Sizes for static/power tilt/power recline seating:


- Support heights: 13", 16", 18", 20", 22"
- Widths: 14", 15", 16", 17", 18", 19", 20", 22", 24"
- Additional back sizes are available with special order





³ OF BACK SUPPORT HEIGHT ADJUSTMENT BUILT INTO THE TRU-BALANCE® 3 BACK SPINE.

TRU-Comfort 2 Back

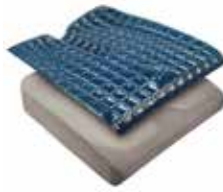

- Adjustable positioning and pressure management
- Compatible with both TRU-Balance® 3 and TRU-Balance® 4 Power Positioning Systems and the static seat frame system
- Four-way stretch, recovery, fluid-resistant, breathable cover



Premiere P	Glacial SP	Zen SP	Simplicity G
 <p>Gray X-Soft SunMate FRG foam topper, high density positioning wedges built in for offloading.</p>	 <p>New SIM G base, fluid resistant overlay, top and sides</p> <p>TPU heat-sealed, low profile, lightweight cooling gel topper</p>	 <p>New SIM G base, SunMate FRG XXSoft (blue topper), low profile, lightweight</p>	 <p>High-density foam</p> <p>Stealth-Tek cover standard</p> <p>More contour to sacral/coccyx relief well</p>
<p>14-15" W: 79.4 kg 16-21" W: 158.8 kg 22" W: 226.8 kg</p>	<p>10-15" W: 79.4 kg 16-21" W: 158.8 kg 22+" W: 181.4 kg</p>	<p>10-15" W: 79.4 kg 16-21" W: 158.8 kg 22+": 181.4 kg</p>	<p>10-15" W: 79.4 kg 16-21" W: 158.8 kg 22+" W: 181.4 kg</p>
<p>1.6 kg</p>	<p>1.6 kg</p>	<p>1.1 kg</p>	<p>.9 kg</p>
<p>8.9 cm / 3.5"</p>	<p>6.4 cm / 2.5"</p>	<p>7.6 cm / 3"</p>	<p>5 cm / 2"</p>
<p>14x13" up to 22x22"</p>	<p>10x10" up to 24x24"</p> <p>Odd sizes in 15-21" W</p>	<p>10x10" up to 24x24"</p> <p>Odd sizes in 15-21" W</p>	<p>10x10" up to 24x24"</p> <p>Odd sizes in 15-21" W</p>

1) CoolCore® is standard on all cushions except the Simplicity G.

UniLink™ Hardware System

- New system from Stealth Products® for mounting thigh guides and elbow blocks onto the TRU-Balance® 3 Seating System
- Elbow blocks in 3"x10", 3"x3", and 4"x4" pad sizes are made with patented Coolcore® technology cover material
- Thigh guides are available in 4"x4", 4"x6", 4"x8", 4"x10", and 4"x12" pad sizes

THIGH GUIDE OPTIONS





Q-LOGIC 3 EX
ENHANCED DISPLAY



Q-LOGIC 3e
JOYSTICK CONTROLLER

- Industry-first, magnetic charging port with adapter can be used with any current XLR charger
- Magnetic wall holder keeps charger plug within reach and off floor
- Integrated puddle light provides greater visibility



Q-LOGIC 3 EX
JOYSTICK CONTROLLER

Q-Logic 3 iAccess
Programmable Switchbox



Q-Logic 3 Stand
Alone Joystick



"MY EDGE GIVES ME THE OPPORTUNITY TO USE ADVANCED TECHNOLOGY THAT GIVES ME SO MUCH INDEPENDENCE."

ZOE HERNANDEZ



Q-LOGIC 3

ADVANCED DRIVE CONTROLS

The Q-Logic 3 Advanced Drive Controls System from Quantum Rehab® is like no other power wheelchair electronics system. Designed from clinician and consumer feedback, Q-Logic 3 is user-friendly and highly customisable.

FEATURES

- Wireless Bluetooth® programming through Windows devices
- Bluetooth access to computer, tablet, and phone functions
- iAccess programmable push button and toggle options for customised seating access

Clinic Mode:

- Quickly set up multiple input devices in client evaluations
- Automatically assigned drive modes to input devices



Q-Logic 3 home screen



Q-Logic 3 drive screen

ENHANCED DISPLAY

The optional Q-Logic 3 EX Enhanced Display serves as a separate screen used with specialty drive controls to display system data. With the Enhanced Display's built-in environmental controls, Q-Logic 3 can operate a wide range of appliances that use infrared (IR) remote controls, such as televisions, and serves as a Bluetooth® computer mouse emulator.



Q-Logic 3 Specialty Controls Interface Module



VR2 6 KEY JOYSTICK



JOYSTICK OPTIONS



SIP-N-PUFF



Stealth® i-Drive®
Sip-N-Puff

Q-LOGIC 3



Q-Logic 3 EX
Joystick



Q-Logic 3e
Joystick



Q-Logic 3 EX
Enhanced Display



Q-Logic 3
Stand Alone
Joystick



Q-Logic 3 iAccess
Programmable
Switchbox



Q-Logic 3
Specialty Interface
Module

PROPORTIONAL DRIVE CONTROLS



Mo-Vis™
All-Round
Joystick



Mo-Vis™
All-Round Light
Joystick



Mo-Vis™
Multi Joystick



Stealth® Precision
Mini Proportional
Joystick



Stealth®
Mushroom
Joystick



Mo-Vis™
Micro Joystick



Mo-Vis™
Multi-Swing 2G



Movis Chin
Control Harness

MULTI-SWITCH

ACCESSORIES



Stealth® i-Drive® 4 Stealth Pro,
3-Switch Head Array



Stealth® i-Drive® Tri-Array,
3-Switch Head Array



Stealth® i-Drive®
Curved Array

SINGLE SWITCH



Stealth®
Egg Switch



Mo-Vis™
Twister Basic
Button



Mo-Vis™
Twister Satellite
Button

TRAY ARRAYS



Stealth® i-Drive®
Eclipse Tray with
Four Capacitive
Proximity Sensors



Stealth® i-Drive®
Eclipse Tray with
Four Fiber Optic
Proximity Sensors

MISC. ACCESSORIES



Front and Rear
LED Fender Lights



Hydration System



Transfer Bars (Pair)



Quantum® Backpack



Personal Item Hook



Phone Holder


































Reversing Camera



Glove Box

ADDITIONAL ACCESSORIES ARE AVAILABLE ON SCRIPT ORDER FORMS. CONSULT YOUR PROVIDER TO DISCUSS/REVIEW OPTIONS.

COLOUR GUIDE

		STRETTO®	EDGE® 3	Q6 EDGE® 2.0	4FRONT® 2	EDGE® 4	Q6 EDGE® HD	Q1450	J623 2.0	R-TRACK®	OUTBACK® 4X4
	Deep Purple	•	•	•	•		•			•	
	Cosmic Purple					•					
	Sugar Plum								•		
	Midnight Blue		•	•	•		•	•		•	
	Blue Jean	•	•			•					
	Electric Blue			•	•		•	•	•	•	
	Steel Blue					•					
	Royale Blue					•					
	Military Green		•								
	Night Train		•	•	•	•					
	Bermuda Blue					•					
	Back in Black	•	•	•	•	•	•	•	•	•	
	Grey Street			•							
	White Satin			•							
	Ghost Gray					•					
	Green Machine	•	•	•			•			•	
	Forest Green					•					
	Lemon Crush		•	•		•				•	
	Evolution Orange		•	•	•	•				•	
	Red Corvette	•	•	•	•		•	•	•	•	
	Inferno Red					•					
	Cherry Bomb		•								
	Raspberry Beret	•	•								
	Pretty in Pink			•		•			•	•	
	Ice Pink					•					
	Lime Light				•				•		
	Sage (Matte)										•
	Tan (Matte)										•
	Camo - Black										•
	Camo - White										•
	Camo - Desert										•

* ADDITIONAL J623 2.0 COLOURS INCLUDE: CARBON FIBRE, EMERALD GREEN, ORANGE CRUSH, SILVER BULLET, & YELLOW BLAZE. CONTACT YOUR NEAREST AUTHORISED QUANTUM PROVIDER TO LEARN MORE.

TRU BALANCE[®] 4

POWER POSITIONING SYSTEMS



ELEVATION



POWER RECLINE



POWER TILT



Available with
12" iLevel



Available with
12" iLevel



Available with
12" iLevel



	EDGE 3 STRETTO®	EDGE® 3	EDGE® 4
DRIVE CONFIGURATION	Mid-Wheel Drive	Mid-Wheel Drive	Mid-Wheel Drive
DRIVE WHEELS	35.5 cm (14") or 31.7 cm (12.5") Pneumatic or Flat Free Tyres	35.5 cm (14") Pneumatic or Flat Free Tyres 14x4" wide Flat Free Tyres ²	35.5 cm (14") Pneumatic or Flat Free Tyres 14x4" wide Flat Free Tyres ²
CASTOR WHEELS	Front: 12.7 cm (5") or 15.2 cm (6") Solid Rear: 15.2 cm (6") Solid	Front: 15.2 cm (6") Solid Rear: 15.2 cm (6") Solid	Front: 17.8 cm (7") Solid Rear: 17.8 cm (7") Solid
ANTI-TIP WHEELS	N/A	N/A	N/A
MAXIMUM SPEED³	Up to 10 km/h (5.6 km/h w/ iLevel) ¹⁰	Up to 10 km/h (7.2 km/h with iLevel)	Up to 10 km/h (7.2 km/h with iLevel)
GROUND CLEARANCE	7.3 cm (2.9")	6.8 cm (2.7") (battery tray)	7.28 cm (2.87") (battery tray)
TURNING RADIUS⁴	53.3 cm (21")	52 cm (20.5")	57.9 cm (20.8")
OVERALL LENGTH⁸	82.5 cm (32.5")	90 cm (35.5") without foot riggings	97.2 cm (22.8") without foot riggings
BASE WIDTH	54.6 cm (21.5") (35.5 cm (14") drive wheels) 51.9 cm (20.47") (31.5 cm (12.5") drive wheels)	61.9 cm (24.4")	61 cm (24.2")
SEAT SIZES	TRU-Balance® 3 and 4	TRU-Balance® 3 and 4	TRU-Balance® 3 and 4
TRU-BALANCE® 4	W: 16-20" D: 16-20"	W: 16-22" D: 16-22"	W: 16-22" D: 16-24"
TRU-BALANCE 3® STATIC	W: 12-20" D: 12-20"	W: 12-22" D: 10-22"	W: 12-22" D: 12-24"
TRU-BALANCE® 3 POWER TILT	W: 12-20" D: 12-20"	W: 12-22" D: 12-22"	W: 12-22" D: 12-24"
TRU-BALANCE® 3 LIFT & TILT	W: 12-20" D: 12-20"	W: 12-22" D: 12-22"	W: 12-22" D: 12-24"
SEAT-TO-FLOOR HEIGHTS^{13,14}	TRU-Balance® 3 and 4	TRU-Balance® 3 and 4	TRU-Balance® 3 and 4
TRU-BALANCE® 4	17.5-19.5"	17.5-19.5"	18.25-21.25"
SYNERGY / STATIC	16-18"	16.5-18.5"	17.75-19.75"
TRU-BALANCE® 3 POWER TILT⁷	17.5-19.5"	17.5-19.5"	17.5-19.5"
TRU-BALANCE® 3 LIFT & TILT⁷	17.5-19.5"	17.5-19.5"	17.5-19.5"
BATTERY SIZE⁵	MK 48Ah(2)	MK 50Ah(2)	MK 70Ah(2)
BATTERY WEIGHT¹¹	17 kg each	17.2 kg each	23.3 kg each
AVAILABLE ELECTRONICS⁹	120A Q-Logic 3	120A Q-Logic 3	120A Q-Logic 3
BATTERY CHARGER	8A Off-Board	8A Off-Board	8A Off-Board
MOTOR PACKAGES	2-Pole	2-Pole or 4-Pole	4-Pole
WEIGHT CAPACITY⁶	136 kg	136 kg	150 kg
APPROX. TOTAL WEIGHT¹⁵	185 kg	193 kg	193 kg

Available with
12° iLevelAvailable with
10° ElevationAvailable with
10° ElevationAvailable with
10° ElevationColour shown
not available

No Elevation



Q6 EDGE® 2.0	4 FRONT® 2	4 FRONT® 2 HD	Q6 EDGE® HD		QUANTUM® 1450
Mid-Wheel Drive	Front-Wheel Drive	Front-Wheel Drive	Mid-Wheel Drive		Front-Wheel Drive
35.5 cm (14") Pneumatic or Flat Free Tyres 14x4" wide Flat Free Tyres ²	35.5 cm (14") Pneumatic or Flat Free Tyres 14x4" wide Flat Free Tyres ²	35.5 cm (14") Pneumatic or Flat Free Tyres 14x4" wide Flat Free Tyres ²	35.5 cm (14") Pneumatic or Flat Free Tyres 14x4" wide Flat Free Tyres ²		35.5 x 10 cm (14x4" wide) Flat Free Tyres
Front: 15.2 cm (6") Solid ¹ 20.3 cm (8") Solid ² Rear: 15.2 cm (6") Solid ¹ 20.3 cm (8") Solid ²	Rear: 22.8 cm (9") Solid	Rear: 22.8 cm (9") Solid	Front: 15.2 cm (6") Solid Rear: 15.2 cm (6") Solid		22.8 cm (9") Solid
N/A	Front: 10 cm (4") Solid	Front: 10.1 cm (4") S	N/A		Front: 7.5 cm (3") Solid
Up to 10 km/h (5.6 km/h w/ iLevel) ¹⁰	Up to 10 km/h (5.6 km/h w/ elevation)	Up to 8 km/h (4.3 km/h w/ elevation)	Up to 7.6 km/h (5.6 km/h w/ iLevel)		Up to 7.2 km/h
6.8 cm (2.7") (battery tray)	7.6 cm (3") (battery tray)	7.6 cm (3") (battery tray)	7.6 cm (3") (frame centre)		5 cm (2") (anti-tips)
52 cm (20.5") w/ 15.2 cm (6") castors 59 cm (23.2") w/ 20.3 cm (8") castors	65.7 cm (25.9")	65.7 cm (25.9")	50 cm (19.7")		67.3 cm (26.5")
90 cm (35.5") w/ 15.2 cm (6") castors 102.6 cm (40.4") w/ 20.3 cm (8") castors	101 cm (39.9")	101 cm (39.9")	90.4 cm (35.6")		97 cm (38.2")
61.9 cm (24.4")	60.9 cm (24")	60.9 cm (24")	66.5 cm (26.2")		74 cm (29.2")
TRU-Balance® 3	TRU-Balance® 3 and 4	TRU-Balance® 3	TB2 HD	TB3 HD	TRU-Balance® 2 HD Tilt
N/A	W: 16-22" D: 16-22"	N/A	N/A	N/A	N/A
W: 12-22" D: 12-22"	W: 12-22" D: 12-22"	HD W: 20-24" D: 20-22"	W: 20-32" D: 16-25"	W: 20-32" D: 16-25"	W: 20-32" D: 20-28"
W: 12-22" D: 12-22"	W: 12-22" D: 12-22"	HD W: 20-24" D: 20-22"	W: 20-32" D: 16-25"	W: 20-32" D: 16-25"	W: 20-32" D: 20-28"
W: 12-22" D: 12-22"	W: 12-22" D: 12-22"	HD W: 20-24" D: 20-22"	N/A	W: 20-32" D: 16-25"	N/A
TRU-Balance® 3	TRU-Balance® 3 and 4	TRU-Balance® 3	TB2 HD	TB3 HD	TRU-Balance® 2 HD Tilt
17.5-19.5"	17.75-19.75"	N/A	N/A	17.7-19.7"	N/A
16.5-18.5"	15.2-17.2"	15-17"	16.5-18.5"	17.5-19.5"	16.5-18.5"
17.5-19.5"	17.75-19.75"	17.75-19.75"	18.2-20.9"	17.7-19.7"	18-20"
17.5-19.5"	17.75-19.75"	17.75-19.75"	N/A	18.5-20.5"	N/A
MK 50Ah(2)	MK 70Ah(2)	MK 70Ah(2)	MK 70Ah(2)		MK 70Ah(2)
17.2 kg each	23.3 kg each	23.3 kg each	23.5 kg each		23.5 kg each
120A Q-Logic 3	120A Q-Logic 3	120A Q-Logic 3	120A Q-Logic 3		90A VR2 6-Key 120A Q-Logic 3
8A Off-Board	8A Off-Board	8A Off-Board	8A Off-Board		8A Off-Board
2-Pole	4-Pole	4-Pole	4-Pole ¹		4-Pole
136 kg	136 kg	204 kg	204 kg ¹²		272 kg
192 kg	201 kg	210 kg	216 kg		203 kg

Available with
10° Elevation



Available with
10° Elevation



Available with
10° Elevation



	J623 2.0	R-TRAK	QUANTUM OUTBACK 4x4
DRIVE CONFIGURATION	Mid-Wheel Drive	Rear-Wheel Drive	4 Wheel Drive
DRIVE WHEELS	35.5 cm (14") Pneumatic or Flat Free Tyres 14x4" wide Flat Free Tyres ²	35.6 cm (14") Pneumatic or Solid Tyres 14x4" wide Flat Free Tyres ²	38.1 cm (14") Knobby Pneumatic Tyres
CASTOR WHEELS	Front: 15.2 cm (6") Solid ¹ 20.3 cm (8") Solid ² Rear: 15.2 cm (6") Solid ¹ 20.3 cm (8") Solid ²	Front: 22.9 cm (9") Solid or Pneumatic	N/A
ANTI-TIP WHEELS	N/A	Rear: 10 cm (4") Solid	N/A
MAXIMUM SPEED³	Up to 10 km/h (2.5 km/h w/ elevation)	Up to 10 km/h (5.6 km/h w/ elevation)	Up to 10 km/h (5.6 km/h w/ elevation)
GROUND CLEARANCE	6.9 cm (2.7") (battery tray)	10 cm (3.97") (battery tray)	10.7 cm (4.2") (motors)
TURNING RADIUS⁴	52 cm (20.5")	65.3 cm (25.7")	125 cm (49.2")
OVERALL LENGTH⁸	90 cm (35.5") without foot riggings	88.3 cm (34.8")	106.7 cm (42")
BASE WIDTH	61.9 cm (24.4")	61.4 cm (24.2")	73.2 cm (28.8")
SEAT SIZES	TRU-Balance® 3	TRU-Balance® 3	TRU-Balance® 3
TRU-BALANCE® 4	N/A	N/A	N/A
TRU-BALANCE 3® STATIC	W: 12- 22" D: 12-22"	W: 12- 24" D: 12-22"	N/A
TRU-BALANCE® 3 POWER TILT	W: 12- 22" D: 12-22"	W: 12- 24" D: 12-22"	W: 16- 24" D: 16-22"
TRU-BALANCE® 3 LIFT & TILT	N/A	W: 12- 24" D: 12-22"	W: 16- 24" D: 16-22"
SEAT-TO-FLOOR HEIGHTS^{13,14}	TRU-Balance® 3	TRU-Balance® 3	TRU-Balance® 3
TRU-BALANCE® 4	17.5-19.5"	17.2-19.2"	N/A
SYNERGY / STATIC	16.5-18.5"	16-18"	N/A
TRU-BALANCE® 3 POWER TILT⁷	17.5-19.5"	17.2-19.2"	19-20"
TRU-BALANCE® 3 LIFT & TILT⁷	17.5-19.5"	17.2-19.2"	20-21"
BATTERY SIZE⁵	MK 50Ah(2)	MK 70Ah(2)	MK 88Ah(2)
BATTERY WEIGHT¹¹	17.2 kg each	23.3 kg each	28 kg each
AVAILABLE ELECTRONICS⁹	90A VR2 6-Key	120A Q-Logic 3	120A Q-Logic 3
BATTERY CHARGER	8A Off-Board	8A Off-Board	8A Off-Board
MOTOR PACKAGES	2-Pole	4-Pole	4-Pole
WEIGHT CAPACITY⁶	136 kg	150 kg	160 kg
APPROX. TOTAL WEIGHT¹⁵	185 kg	201 kg	226 kg

1 Standard. **2** Optional. **3** Speed varies with user weight, terrain type, battery charge, battery condition and tyre pressure. **4** Without front riggings. **5** Battery capacity may change when adding power positioning systems. **6** Weight capacity is dictated by the selection of seating. **7** All measurements +/- 1/8" with 18" x 18" seat, fixed towers, 136 kg weight capacity and solid tyres (add 3/8" with pneumatic tyres). Depths may increase seat-to-floor height. Contact Customer Service for details. Listed measurements are in their lowest position; does not include elevating system unless specified. **8** Without foot riggings. **9** All Quantum units can be ordered with specialty controls. **10** 10 km/h speed availability only with Q-Logic 3. **11** Battery weight may vary +/- 1.3 kg based on manufacturer. **12** Weight capacity is reduced to 181 kg with iLevel® power adjustable seat height. **13** Seat-to-floor heights are measured with the unit unloaded. Due to the 4Front's suspension system, seat-to-floor heights are approximately 0.6 cm lower when the unit is loaded. **14** 41.9 cm seat-to-floor height is not available with power tilt or iLevel®. 49.5 cm seat-to-floor height is not available with static seat or power recline. **15** these values are with batteries installed, TB3 Tilt, Lift, Recline and Centre Mount Power Elevating Leg Rests. Overall maximum weight will change depending on Seating scripted. **NOTE:** The information contained herein is



GET AMONGST IT! with Stocky

WATCH THE COMPLETE SERIES ON



PQ365 APP


ALL QUANTUM® PRODUCT INFORMATION IN ONE APP!

Designed for providers, the free, easy-to-use app helps busy professionals stay up to date with the latest information and allows access all relevant content in one quick and convenient place.

The app has many cool features, including, serial and order number search, order delivery status, invoice number look-up, a voice search function, an image share option, technical support, marketing material and educational resources. Additional features include a quick-search option, 360-degree product views, push notifications, news, events, specifications sheets, brochures, owner's manuals, basic operating instructions and more.



Download PQ365 today!

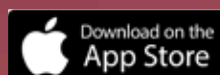



QUANTUMREHAB.COM.AU

The official Australian website for Quantum® is fully equipped with the latest Australian product information and resources along with valuable provider links.



RESOURCEHUB.PRIDEMOBILITY.COM.AU

The Quantum Resource Hub is a platform designed to improve your experience and streamline your business operations with us. It offers a variety of resources to help manage your business more efficiently and effectively.

QUICK TIPS WITH QUANTUM

Your go to for power chair tips, tricks and clinical hacks that make your day smoother and scripting smarter! Watch on our social channels.



QUANTUMREHAB.COM

Quantum's official website is the ultimate source for everything you need to know about Quantum® products. Highlights include detailed product specifications, order forms, videos, user testimonials, and much more.



Q3DC.QUANTUMREHAB.COM (ECON-W)


Quantum Rehab's® login page for the Q-Logic 3 cloud-based wireless programming. Q-Logic 3 is our newest generation of electronics available with improved features and benefits inspired by consumer, provider and therapist feedback.



SOCIAL MEDIA

Stay up to date on the latest Quantum® events and product releases through LinkedIn, YouTube, Facebook, and Instagram.

SERVICE AND SUPPORT

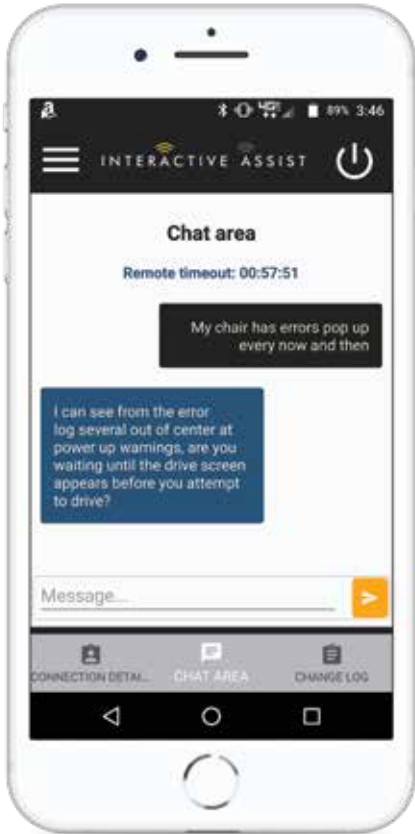

INTERACTIVE ASSIST

Q-Logic 3 now features Interactive Assist. Interactive Assist from Quantum Rehab® is an app which provides a direct real-time connection from a power chair's electronics to a provider's programming station. The connection provides the technician with full system and diagnostics information along with tools, such as a real-time mirror image of the power chair's electronics display, to eliminate wasteful trips.


- Change Log History
- Saving XPA Files
- Motor Usage History
- Text Chatting
- Programming IR Menu
- Access to Bluetooth Features
- Dealer-to-Dealer Programming Access
- Screen Mirroring
- Real-time Monitoring
- Efficiency Improvements
- Troubleshooting
- Battery Management
- Diagnostic History
- Basic ProgramminG
- Remote Access



Download Interactive Assist today!

ECON - I

The Econ-I app is now available for download on Apple® iPhone® devices, allowing you to program Quantum® Power Chairs conveniently from your iPhone or iPad. Compatible devices include the iPhone 6S/6s Plus/SE or higher, and the iPad Air/iPad Mini 2/iPad Pro/iPad 5th generation or higher. All devices must support iOS version 12 or higher.

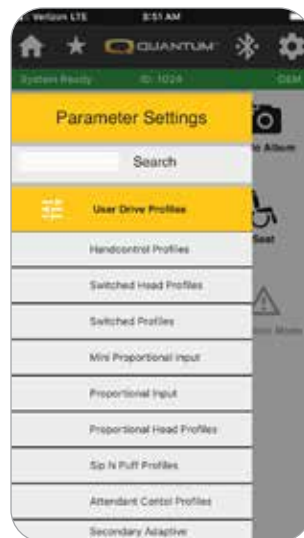

Download ECON - I today!



EPiC seating

effortless postural control



Back and Head Support sold separately. But any after-market back that can be mounted with a 1in (2.54cm) cane clamp hardware will fit your new **EPiC Seating**.

Seating setup as shown:

Head Support:
Comfort Plus™

Head Support Hardware:
UniLink™

Back:
ADI® Aluminum Back, also check TFB™ Contour with similar configurations.

Back Mounting:
Stealth Quick Release/Pro
Simple Solutions™ FlipLock™

Fits Frame Width: in(cm)

Small	10 to 16 (25.40 to 40.64)
Medium	16 to 21 (40.64 to 53.34)
Large	21 to 26 (53.34 to 66.04)

check EPiC Seating compatible frames

1in (2.54cm)
Cane Diameter

Quick Adjustment Lever
adjust anatomical pivot-point leaning angle without the need of hardware

Load Absorbers
come in two resistance levels

50N ideal for users with limited active movement or posterior pelvic tilt

100N ideal for users with normal to highly active movement

Initial Angle Setup
personalized initial seating frame setup points

Pelvic Pivot-Point Mount
rotates around an imaginary radius to the user's acetabulum or the greater trochanter.



6in (15.24cm)
posterior pelvic support height*


* plate height, more options coming soon



Posterior Pelvic Support



EPiC Seating Pivot Reference

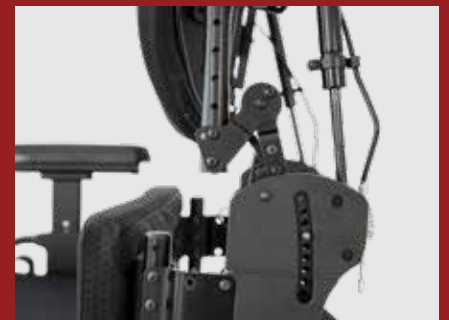

Pivot Mechanics

STEALTH RESOURCES

STEALTHPRODUCTS.COM



Stealth Products® offers a range of versatile products, including head and body support/positioning, joystick mounts, trays, seating, cushions, and paediatric mobility products.



CUSHIONS.STEALTHPRODUCTS.COM/

Learn about our Stealth Products range of cushions featuring our exclusive Coolcore® Technology.



Scan here to access the Stealth Cushions Catalogue



ADI.STEALTHPRODUCTS.COM/

Learn about our Stealth Products range of ergonomically designed ADI® back supports and seating systems.



Scan here to access the ADI Catalogue



UNILINK.COM

Learn about our Stealth Products UniLink range of strong, unique, and simplified components.



Scan here to access the UniLink Catalogue



HEADSOLUTIONS.COM

Learn about our Stealth Products Head Solutions range of versatile, high quality head and neck positioning systems and mount options.



Scan here to access the Head Solutions Catalogue

quantumrehab.com.au



QUANTUM®

1 FOR REHAB POWER

